




Waterbath

WNB 7

The combination of corrosion-resistant stainless steel, precise electronics and multiple temperature protection guarantees the highest level of safety in the laboratory.



On this page, you can find all the essential technical data on the Memmert water bath. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at myAtmoSAFE@memmert.com.

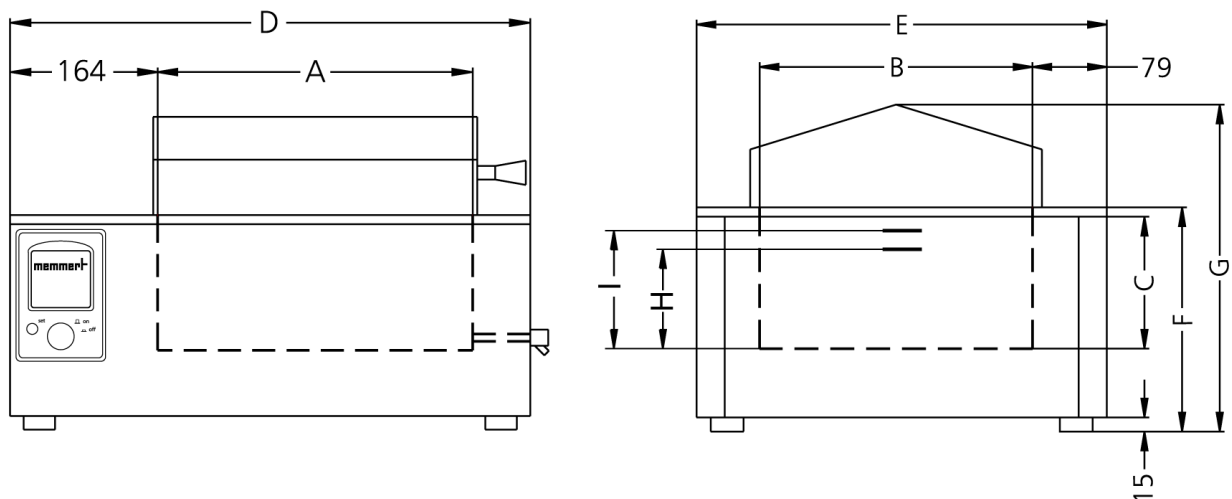

Temperature

Temperature range in °C	min. 5°C above ambient up to +95°C with additional boiling mode (+100°C)
-------------------------	--

Control of standard components

Temperature	1 Pt100 sensor class A in 4-wire-circuit
Timer	integrated digital timer from 1 min. to 99,99 hours for:ON continuous operationWAIT (delayed on for continuous andlimited timed operation)HOLD

Safety

Temperature control	mechanical temperature limiter TB protectionclass 1 switching the heating off at approx. 30°Cabove max. temperature of the bath
Temperature control	in case of overtemperature due to failure,the heating is switched off at approx. 10°Cabove the set temperature (fixed value)
Autodiagnostic system	microprocessor PID-temperature controllerwith integrated autodiagnostic system withfault indicator

Heating concept

Heating Baths	corrosion-proof large-area heating onthree sides
---------------	--

Standard equipment

Controller	digital display (LED) of set and actualtemperature (0,1°C resolution) and of (remaining)programme time
------------	--

Stainless steel interior

Dimensions W x H x D in mm	$w_{(A)} \times h_{(C)} \times d_{(B)}$: 240 x 140 x 210 mm
Interior	easy-to-clean interior, made of stainlesssteel, reinforced by deep drawn ribbing,material no. 1.4301 (ASTM 304), laser-welded
Volume	7 l
Liquid level min.	(H) 97 mm
Liquid level max.	(I) 120 mm

Textured stainless steel casing

Dimensions	$w_{(D)} \times h_{(G)} \times d_{(E)}$: 468 x 356 x 238 mm
------------	--

Electrical data

Voltage	230 V, 50/60 Hz
Electrical load	approx. 1200 W (during heating)

Ambient conditions

Installation	The vent openings in the left and back side must remain unobstructed. Minimum wall spacing on all sides is 80 mm. The minimum spacing from the top of the bath to the next ceiling is 750 mm.
Ambient temperature	+5 °C to +40 °C
Humidity rh	max. 80 %, non-condensing
Overvoltage category	II
Pollution degree	2

Packing/shipping data

Transport information	The appliances must be transported upright
Customs tariff number	8419 8998
Country of origin	Federal Republic of Germany
WEEE-Reg.-No.	DE 66812464
Dimensions approx incl. carton	B x H x T: 570 x 450 x 420 mm
Net weight	approx. 11 kg
Gross weight carton	approx. 16 kg

Standard units are safety-approved and bear the test marks

